

BORIS

The comfortable toilet and shower chair



BORIS – the comfortable toilet and shower chair

- The anatomically moulded, deep seat pan made of PE plastic with torso and side guides provides particularly high seat stability.
- BORIS offers a robust and easy-care subframe made of stainless steel.
- Height adjustable seat pan for adaptation to the respective constructional situation.
- The seat angle is steplessly adjustable via a gas pressure spring from (-)5° to 30°.
- For support in sitting position, the BORIS toilet and shower chair can be combined with our fixation system FLEXI-AKTIV.
- We offer head rests and foot supports, insert and arm support pads and a plug-in spray protection as accessories for individual adaptation.
- A larger hygiene opening is available for sizes 2 - 4.



Anatomically moulded seat pan with extremely good side guides for stabilisation of the upper body. Two-part insert pad for additional sitting comfort



Height and angle adjustable head rest for support of head control



Foldable and angle adjustable foot support made of high-quality ABS material




Plug-in arm support for stabilisation and for keeping the upper body upright




Plug-in spray protection available in two sizes



Practical toilet arrangement

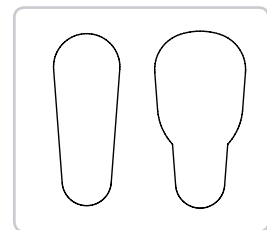

Larger hygiene opening available for BORIS size 2 - 4

	Size 0	Size 1	Size 2	Size 3	Size 4
 Seat width at front	31 cm	36 cm	41 cm	46 cm	51 cm
Seat width at rear	26 cm	31 cm	36 cm	41 cm	46 cm
Seat depth	31 cm	36 cm	41 cm	46 cm	48 cm
Back height	49 cm	51 cm	56 cm	61 cm	65 cm
Seat angle	(-)5° - 30°	(-)5° - 30°	(-)5° - 30°	(-)5° - 30°	(-)5° - 30°
Can be positioned over a toilet	47 - 57cm	45 - 55cm	42 - 52cm	42 - 52cm	42 - 52cm
Total width	56 cm	61 cm	63 cm	71.5 cm	74 cm
Total length	65 cm	65 cm	65 cm	70 cm	70 cm
Max. lower leg length	38 cm	38 cm	48 cm	48 cm	48 cm
Front edge of hygiene opening to back edge of subframe	47 cm	49 cm	51 cm	54 cm	57 cm
Width hygiene opening at front	10 cm	10 cm	10 cm	10 cm	10 cm
Width hygiene opening at back	12 cm	12 cm	13.5 cm	13.5 cm	13.5 cm
Width hygiene opening at back*	-	-	18.5 cm	18.5 cm	18.5 cm
Max. load	60 kg	80 kg	100 kg	100 kg	150 kg
HMV No.	33.40.04.2008				



* only applies to the sizes 2 b, 3 b and 4 b

Frame colour:  Stainless steel

BORIS